

## NB1-J02 - October 13, 2015

Item # NB1-J02 was discontinued on October 13, 2015. For informational purposes, this is a copy of the website content at that time and is valid only for the stated product.

### ARGON ION LASER MIRRORS

- ▶ 1 kW/cm<sup>2</sup> Damage Threshold
- ▶ Designed for Either 0° or 45° AOI
- ▶ R<sub>avg</sub> > 99.5% Over Wavelength Ranges



**NB1-H07**

Ø1" Ar-Ion Laser Line Mirror, 0° AOI



**NB1-J05**

Ø1" Ar-Ion Laser Line Mirror, 45° AOI



**NB07-J02**

Ø19 mm Ar-Ion Laser Line Mirror, 45° AOI



Back Side  
Engraved with Item #

[Hide Overview](#)

#### OVERVIEW

##### Features

- Available Wavelength Ranges
  - 244 - 257 nm (45° AOI)
  - 300 - 308 nm (45° AOI)
  - 333 - 364 nm (0° and 45° AOI)
  - 458 - 528 nm (45° AOI)
- Ø19 mm and Ø1" Options

Ar-Ion laser line mirrors have a specialized coating that provides a high reflectance over wavelength ranges coincident with the emission of an Ar-Ion laser. The coating has a high damage threshold suitable for use with the output beam from large-frame high-powered CW (continuous wave) Ar-Ion lasers.

Our selection of Ø19 mm mirrors are specifically designed to fit our Polaris Fixed Optic Mounts for laser system design and other OEM applications. This diameter provides a larger clear aperture than Ø1/2" optics while allowing the mounts to maintain a Ø1" footprint.

Specification	Value
Average Reflectance	>99.5%
Damage Threshold at 45° AOI	1 kW/cm <sup>2</sup>
Substrate Material	Fused Silica
Front Surface Flatness	λ/10
Clear Aperture	>80% of Diameter
Back Surface	Fine Ground <sup>a</sup>
Parallelism	≤3 arcmin
Diameter Tolerance	+0.0/-0.1 mm
Thickness	6.0 ± 0.2 mm
Surface Quality	10-5 Scratch-Dig

a. A fine ground back surface is frosted and will diffuse light that is not reflected by the mirror's front surface.



[Hide Specs](#)

## S P E C S

General Specifications	
Substrate Material	Fused Silica
Front Surface Flatness	$\lambda/10$
Clear Aperture	>80% of Diameter
Back Surface	Fine Ground <sup>a</sup>
Diameter Tolerance	+0.0/-0.1 mm
Thickness	6.0 ± 0.2 mm
Surface Quality	10-5 Scratch-Dig
CW Damage Threshold	1 kW/cm <sup>2</sup>

a. A fine ground back surface is frosted and will diffuse light that is not reflected by the mirror's front surface.

Ar-Ion ( $\lambda = 244-257$  nm)

Item #	NB07-J02	NB1-J02
Diameter	Ø19.0 mm	Ø25.4 mm
Reflectance	$R_{avg} >99.5\%$	
Angle of Incidence	45°	

Ar-Ion ( $\lambda = 300-308$  nm)

Item #	NB07-J05	NB1-J05
Diameter	Ø19.0 mm	Ø25.4 mm
Reflectance	$R_{avg} >99.5\%$	
Angle of Incidence	45°	

Ar-Ion ( $\lambda = 333-364$  nm)

Item #	NB1-H07	NB07-J07	NB1-J07
Diameter	Ø25.4 mm	Ø19.0 mm	Ø25.4 mm
Reflectance	$R_{avg} >99.5\%$		
Angle of Incidence	0°	45°	45°

Ar-Ion ( $\lambda = 458-528$  nm)

Item #	NB07-J10	NB1-J10
Diameter	Ø19.0 mm	Ø25.4 mm
Reflectance	$R_{avg} >99.5\%$	
Angle of Incidence	45°	

[Hide Damage Thresholds](#)

## D A M A G E T H R E S H O L D S

## Damage Threshold Data for Thorlabs' Argon Ion Laser Mirrors

The specifications to the right are measured data for Thorlabs' argon ion laser mirrors. Damage threshold specifications are constant for a given optical density coating.

Damage Threshold Specifications		
Optical Density Designation (Item # Suffix)	Wavelength Range	Damage Threshold (CW)
-02	244 - 257 nm	1 kW/cm <sup>2</sup>
-05	300 - 308 nm	1 kW/cm <sup>2</sup>
-07	333 - 364 nm	1 kW/cm <sup>2</sup>
-10	458 - 528 nm	1 kW/cm <sup>2</sup>

## Laser Induced Damage Threshold Tutorial

The following is a general overview of how laser induced damage thresholds are measured and how the values may be utilized in determining the appropriateness of an optic for a given application. When choosing optics, it is important to understand the Laser Induced Damage Threshold (LIDT) of the optics being used. The LIDT for an optic greatly depends on the type of laser you are using. Continuous wave (CW) lasers typically cause damage from

thermal effects (absorption either in the coating or in the substrate). Pulsed lasers, on the other hand, often strip electrons from the lattice structure of an optic before causing thermal damage. Note that the guideline presented here assumes room temperature operation and optics in new condition (i.e., within scratch-dig spec, surface free of contamination, etc.). Because dust or other particles on the surface of an optic can cause damage at lower thresholds, we recommend keeping surfaces clean and free of debris. For more information on cleaning optics, please see our *Optics Cleaning* tutorial.


## Testing Method

Thorlabs' LIDT testing is done in compliance with ISO/DIS11254 specifications. A standard 1-on-1 testing regime is performed to test the damage threshold.

First, a low-power/energy beam is directed to the optic under test. The optic is exposed in 10 locations to this laser beam for a set duration of time (CW) or number of pulses (pulse repetition frequency specified). After exposure, the optic is examined by a microscope (~100X magnification) for any visible damage. The number of locations that are damaged at a particular power/energy level is recorded. Next, the power/energy is either increased or decreased and the optic is exposed at 10 new locations. This process is repeated until damage is observed. The damage threshold is then assigned to be the highest power/energy that the optic can withstand without causing damage. A histogram such as that below represents the testing of one BB1-E02 mirror.



Example Exposure Histogram

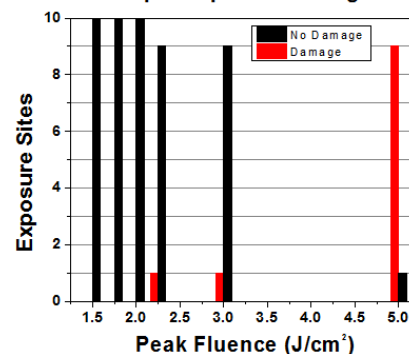

The photograph above is a protected aluminum-coated mirror after LIDT testing. In this particular test, it handled  $0.43 \text{ J/cm}^2$  (1064 nm, 10 ns pulse, 10 Hz,  $\text{Ø}1.000 \text{ mm}$ ) before damage.

According to the test, the damage threshold of the mirror was  $2.00 \text{ J/cm}^2$  (532 nm, 10 ns pulse, 10 Hz,  $\text{Ø}0.803 \text{ mm}$ ). Please keep in mind that these tests are performed on clean optics, as dirt and contamination can significantly lower the damage threshold of a component. While the test results are only representative of one coating run, Thorlabs specifies damage threshold values that account for coating variances.

## Continuous Wave and Long-Pulse Lasers


When an optic is damaged by a continuous wave (CW) laser, it is usually due to the melting of the surface as a result of absorbing the laser's energy or damage to the optical coating (antireflection) [1]. Pulsed lasers with pulse lengths longer than  $1 \mu\text{s}$  can be treated as CW lasers for LIDT discussions. Additionally, when pulse lengths are between 1 ns and  $1 \mu\text{s}$ , LIDT can occur either because of absorption or a dielectric breakdown (must check both CW and pulsed LIDT). Absorption is either due to an intrinsic property of the optic or due to surface irregularities; thus LIDT values are only valid for optics meeting or exceeding the surface quality specifications given by a manufacturer. While many optics can handle high power CW lasers, cemented (e.g., achromatic doublets) or highly absorptive (e.g., ND filters) optics tend to have lower CW damage thresholds. These lower thresholds are due to absorption or scattering in the cement or metal coating.

Pulsed lasers with high pulse repetition frequencies (PRF) may behave similarly to CW beams. Unfortunately, this is highly dependent on factors such as absorption and thermal diffusivity, so there is no reliable method for determining when a high PRF laser will damage an optic due to thermal effects. For beams with a large PRF both the average and peak powers must be compared to the equivalent CW power. Additionally, for highly transparent materials, there is little to no drop in the LIDT with increasing PRF.

In order to use the specified CW damage threshold of an optic, it is necessary to know the following:


1. Wavelength of your laser
2. Linear power density of your beam (total power divided by  $1/e^2$  spot size)
3. Beam diameter of your beam ( $1/e^2$ )
4. Approximate intensity profile of your beam (e.g., Gaussian)

Example Test Data			
Fluence	# of Tested Locations	Locations with Damage	Locations Without Damage
$1.50 \text{ J/cm}^2$	10	0	10
$1.75 \text{ J/cm}^2$	10	0	10
$2.00 \text{ J/cm}^2$	10	0	10
$2.25 \text{ J/cm}^2$	10	1	9
$3.00 \text{ J/cm}^2$	10	1	9
$5.00 \text{ J/cm}^2$	10	9	1

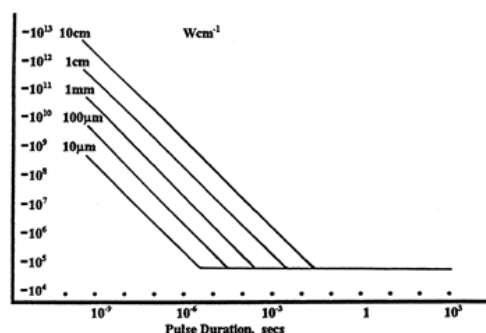

LIDT in linear power density vs. pulse length and spot size. For long pulses to CW, linear power density becomes a constant with spot size. This graph was obtained from [1].

The power density of your beam should be calculated in terms of W/cm. The graph to the right shows why the linear power density provides the best metric for long pulse and CW sources. Under these conditions, linear power density scales independently of spot size; one does not need to compute an adjusted LIDT to adjust for changes in spot size. This calculation assumes a uniform beam intensity profile. You must now consider hotspots in the beam or other nonuniform intensity profiles and roughly calculate a maximum power density. For reference, a Gaussian beam typically has a maximum power density that is twice that of the uniform beam (see lower right).

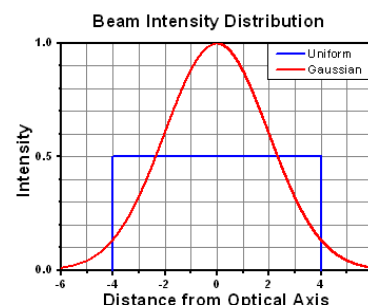

Now compare the maximum power density to that which is specified as the LIDT for the optic. If the optic was tested at a wavelength other than your operating wavelength, the damage threshold must be scaled appropriately. A good rule of thumb is that the damage threshold has a linear relationship with wavelength such that as you move to shorter wavelengths, the damage threshold decreases (i.e., a LIDT of 10 W/cm at 1310 nm scales to 5 W/cm at 655 nm):

$$\text{Adjusted LIDT} = \text{LIDT Power} \left( \frac{\text{Your Wavelength}}{\text{LIDT Wavelength}} \right)$$

While this rule of thumb provides a general trend, it is not a quantitative analysis of LIDT vs wavelength. In CW applications, for instance, damage scales more strongly with absorption in the coating and substrate, which does not necessarily scale well with wavelength. While the above procedure provides a good rule of thumb for LIDT values, please contact Tech Support if your wavelength is different from the specified LIDT wavelength. If your power density is less than the adjusted LIDT of the optic, then the optic should work for your application.

Please note that we have a buffer built in between the specified damage thresholds online and the tests which we have done, which accommodates variation between batches. Upon request, we can provide individual test information and a testing certificate. The damage analysis will be carried out on a similar optic (customer's optic will not be damaged). Testing may result in additional costs or lead times. Contact Tech Support for more information.

### Pulsed Lasers

As previously stated, pulsed lasers typically induce a different type of damage to the optic than CW lasers. Pulsed lasers often do not heat the optic enough to damage it; instead, pulsed lasers produce strong electric fields capable of inducing dielectric breakdown in the material. Unfortunately, it can be very difficult to compare the LIDT specification of an optic to your laser. There are multiple regimes in which a pulsed laser can damage an optic and this is based on the laser's pulse length. The highlighted columns in the table below outline the relevant pulse lengths for our specified LIDT values.


Pulses shorter than  $10^{-9}$  s cannot be compared to our specified LIDT values with much reliability. In this ultra-short-pulse regime various mechanics, such as multiphoton-avalanche ionization, take over as the predominate damage mechanism [2]. In contrast, pulses between  $10^{-7}$  s and  $10^{-4}$  s may cause damage to an optic either because of dielectric breakdown or thermal effects. This means that both CW and pulsed damage thresholds must be compared to the laser beam to determine whether the optic is suitable for your application.

Pulse Duration	$t < 10^{-9}$ s	$10^{-9} < t < 10^{-7}$ s	$10^{-7} < t < 10^{-4}$ s	$t > 10^{-4}$ s
Damage Mechanism	Avalanche Ionization	Dielectric Breakdown	Dielectric Breakdown or Thermal	Thermal
Relevant Damage Specification	N/A	Pulsed	Pulsed and CW	CW

When comparing an LIDT specified for a pulsed laser to your laser, it is essential to know the following:

1. Wavelength of your laser
2. Energy density of your beam (total energy divided by  $1/e^2$  area)
3. Pulse length of your laser
4. Pulse repetition frequency (prf) of your laser
5. Beam diameter of your laser ( $1/e^2$ )
6. Approximate intensity profile of your beam (e.g., Gaussian)

The energy density of your beam should be calculated in terms of J/cm<sup>2</sup>. The graph to the right shows why the energy density provides the best metric for short pulse sources. Under these conditions, energy density scales independently of spot size, one does not need to compute an adjusted LIDT to adjust for changes in spot size. This calculation assumes a uniform beam intensity profile. You must now adjust this energy density to account for hotspots or other nonuniform intensity profiles and

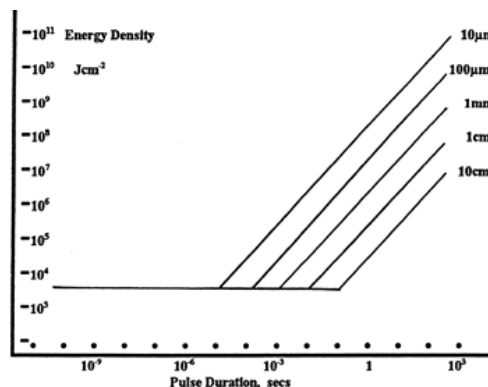

roughly calculate a maximum energy density. For reference a Gaussian beam typically has a maximum energy density that is twice that of the 1/e<sup>2</sup> beam.

Now compare the maximum energy density to that which is specified as the LIDT for the optic. If the optic was tested at a wavelength other than your operating wavelength, the damage threshold must be scaled appropriately [3]. A good rule of thumb is that the damage threshold has an inverse square root relationship with wavelength such that as you move to shorter wavelengths, the damage threshold decreases (i.e., a LIDT of 1 J/cm<sup>2</sup> at 1064 nm scales to 0.7 J/cm<sup>2</sup> at 532 nm):

LIDT in energy density vs. pulse length and spot size. For short pulses, energy density becomes a constant with spot size. This graph was obtained from [1].

$$\text{Adjusted LIDT} = \text{LIDT Energy} \sqrt{\frac{\text{Your Wavelength}}{\text{LIDT Wavelength}}}$$

You now have a wavelength-adjusted energy density, which you will use in the following step.

Beam diameter is also important to know when comparing damage thresholds. While the LIDT, when expressed in units of J/cm<sup>2</sup>, scales independently of spot size; large beam sizes are more likely to illuminate a larger number of defects which can lead to greater variances in the LIDT [4]. For data presented here, a <1 mm beam size was used to measure the LIDT. For beams sizes greater than 5 mm, the LIDT (J/cm<sup>2</sup>) will not scale independently of beam diameter due to the larger size beam exposing more defects.

The pulse length must now be compensated for. The longer the pulse duration, the more energy the optic can handle. For pulse widths between 1 - 100 ns, an approximation is as follows:

$$\text{Adjusted LIDT} = \text{LIDT Energy} \sqrt{\frac{\text{Your Pulse Length}}{\text{LIDT Pulse Length}}}$$

Use this formula to calculate the Adjusted LIDT for an optic based on your pulse length. If your maximum energy density is less than this adjusted LIDT maximum energy density, then the optic should be suitable for your application. Keep in mind that this calculation is only used for pulses between 10<sup>-9</sup> s and 10<sup>-7</sup> s. For pulses between 10<sup>-7</sup> s and 10<sup>-4</sup> s, the CW LIDT must also be checked before deeming the optic appropriate for your application.

Please note that we have a buffer built in between the specified damage thresholds online and the tests which we have done, which accommodates variation between batches. Upon request, we can provide individual test information and a testing certificate. Contact Tech Support for more information.

- [1] R. M. Wood, Optics and Laser Tech. **29**, 517 (1997).
- [2] Roger M. Wood, *Laser-Induced Damage of Optical Materials* (Institute of Physics Publishing, Philadelphia, PA, 2003).
- [3] C. W. Carr *et al.*, Phys. Rev. Lett. **91**, 127402 (2003).
- [4] N. Bloembergen, Appl. Opt. **12**, 661 (1973).


[Hide Ar-Ion Laser Line Mirrors: 244-257 nm](#)

### Ar-Ion Laser Line Mirrors: 244-257 nm



#### Specifications at 45°

- ▶ S-Polarization: R<sub>avg</sub> > 99.5% (244 - 257 nm)
- ▶ P-Polarization: R<sub>avg</sub> > 99.5% (244 - 257 nm)
- ▶ Damage Threshold: 1 kW/cm<sup>2</sup> @ 244 - 257 nm

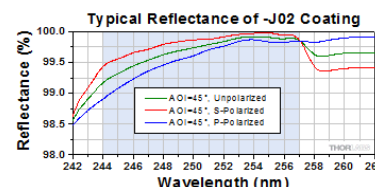

Click to Enlarge  
Click Here for Raw Data

Part Number	Description	Price	Availability
NB07-J02	Ø19 mm Ar-Ion Laser Mirror, 244-257 nm, 45° AOI	\$114.00	Today
NB1-J02	Ø1" Ar-Ion Laser Mirror, 244-257 nm, 45° AOI	\$130.00	Lead Time


[Hide Ar-Ion Laser Line Mirrors: 300-308 nm](#)

### Ar-Ion Laser Line Mirrors: 300-308 nm



#### Specifications at 45°

- ▶ S-Polarization:  $R_{avg} > 99.5\%$  (300 - 308 nm)
- ▶ P-Polarization:  $R_{avg} > 99.5\%$  (300 - 308 nm)
- ▶ Damage Threshold: 1 kW/cm<sup>2</sup> @ 300 - 308 nm

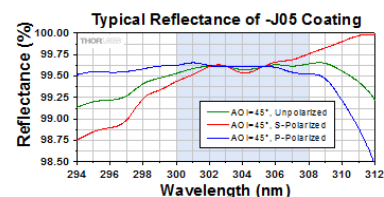

Click to Enlarge  
Click Here for Raw Data

Part Number	Description	Price	Availability
NB07-J05	Ø19 mm Ar-Ion Laser Mirror, 300-308 nm, 45° AOI	\$114.00	Today
NB1-J05	Ø1" Ar-Ion Laser Mirror, 300-308 nm, 45° AOI	\$130.00	3-5 Days

[Hide Ar-Ion Laser Line Mirrors: 333-364 nm](#)

### Ar-Ion Laser Line Mirrors: 333-364 nm




#### NB1-H07 Specifications at 0°

- ▶ Unpolarized:  $R_{avg} > 99.5\%$  (333 - 364 nm)
- ▶ Damage Threshold: 1 kW/cm<sup>2</sup> @ 333 - 364 nm


#### NB1-J07 Specifications at 45°

- ▶ S-Polarization:  $R_{avg} > 99.5\%$  (333 - 364 nm)
- ▶ P-Polarization:  $R_{avg} > 99.5\%$  (333 - 364 nm)
- ▶ Damage Threshold: 1 kW/cm<sup>2</sup> @ 333 - 364 nm

Please note that reflectance does not depend on polarization as the angle of incidence (AOI) approaches 0°.

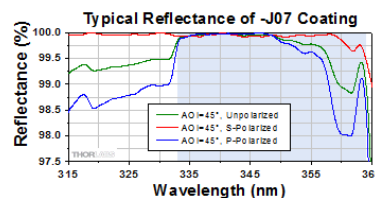


Click to Enlarge  
Click Here for Raw Data

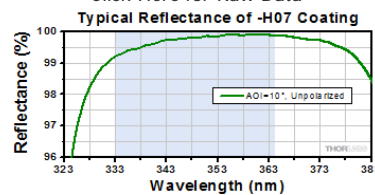

Click to Enlarge  
Click Here for Raw Data


Part Number	Description	Price	Availability
NB1-H07	Ø1" Ar-Ion Laser Mirror, 333-364 nm, 0° AOI	\$126.00	Today
NB07-J07	Ø19 mm Ar-Ion Laser Mirror, 333-364 nm, 45° AOI	\$114.00	Today
NB1-J07	Ø1" Ar-Ion Laser Mirror, 333-364 nm, 45° AOI	\$130.00	Today

## Ar-Ion Laser Line Mirrors: 458-528 nm

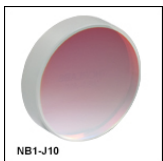

### Specifications at 45°

- ▶ S-Polarization:  $R_{avg} > 99.5\%$  (458 - 528 nm)
- ▶ P-Polarization:  $R_{avg} > 99.5\%$  (458 - 528 nm)
- ▶ Damage Threshold: 1 kW/cm<sup>2</sup> @ 458 - 528 nm

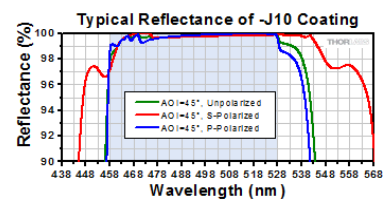

Click to Enlarge  
[Click Here for Raw Data](#)

Part Number	Description	Price	Availability
NB07-J10	Ø19 mm Ar-Ion Laser Mirror, 458-528 nm, 45° AOI	\$114.00	Today
NB1-J10	Ø1" Ar-Ion Laser Mirror, 458-528 nm, 45° AOI	\$130.00	Today