

Indexing Mounting Base



This Indexing Mounting Base is ideal for redirecting a laser to multiple targets on your optical table. Applications include precise switching of shared laser systems, directing optical signals to multiple test instruments, and easy exchange of optical components (with extra top mounts).


For best performance, ensure that the NX1F is securely attached to its support structure. A series of three #8-32 (M4-0.7) mounting holes are provided for non-post mounting applications.



- Rotates in 15 Discrete Steps to Allow Precise Re-Redirecting of Shared Laser
- Repeatability Better than 10 μ rad
- #8-32 (M4) Mounting Holes

ITEM#	METRIC ITEM#	\$	£	€	RMB	DESCRIPTION
NX1F	NX1F/M	\$ 79.00	£ 49.80	€ 73,50	¥ 754.50	Complete 15-Position Indexing Mounting Base
NX1B	NX1B/M	\$ 42.00	£ 26.50	€ 39,10	¥ 401.10	Extra Bottom Plate
NX1FT	NX1FT/M	\$ 37.00	£ 23.30	€ 34,40	¥ 353.40	Extra Top Plate


Kinematic Breadboard Seats

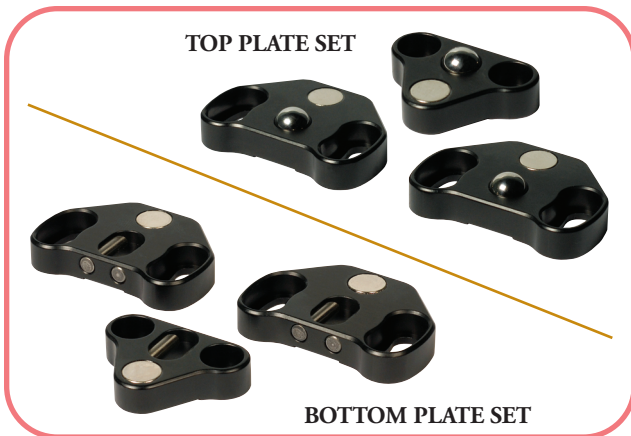

The KBS98 kinematic seats can be used to convert any breadboard into a kinematic platform that can be removed and replaced with a high degree of accuracy. Two of the three seats of each set (top and bottom) have elongated slots to allow the seats to all point to a common center in the middle of the breadboard.

- Three Radial V-Groove Kinematic Design
- Build Your Own Kinematic Base System With Our Breadboards
- Magnetically Coupled With High Strength Rare-Earth Magnets
- 1/4"-20 (M6) 1" (25mm) Centers Turn Any Breadboard Into a Kinematic Base

Slotted design allows the 3 V-Grooves to be aligned to a common center for greatest stability.

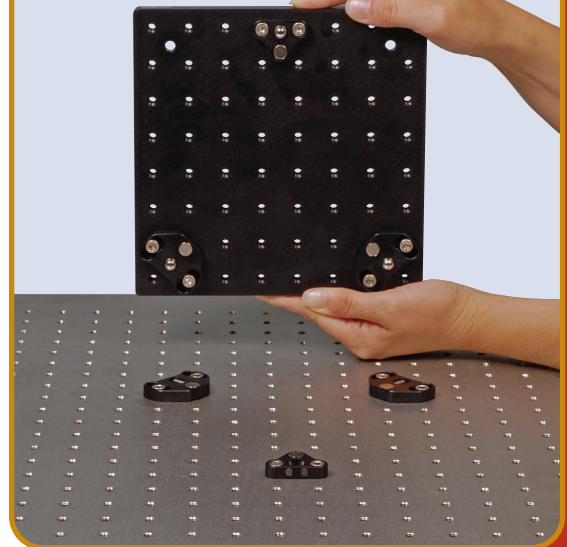

ITEM#	\$	£	€	RMB	DESCRIPTION
KBS98	\$ 133.10	£ 83.90	€ 123,80	¥ 1,271.10	Complete 6 Piece Kinematic Seat Set
KBS98T	\$ 66.55	£ 41.95	€ 61,90	¥ 635.55	3 Piece Top Only Set of Kinematic Seats
KBS98B	\$ 66.55	£ 41.95	€ 61,90	¥ 635.55	3 Piece Bottom Only Set of Kinematic Seats

*Universal Design is Imperial and Metric compatible