

SM3 Endcap

The SM3CP2 endcap can be used to protect delicate optical assemblies from dust when not in active use. The center of the endcap has a drill point to facilitate custom machining.



SM3CP2

Also Available
in Ø1/2",
Ø1.00", Ø2.00"

ITEM#	\$	£	€	RMB	DESCRIPTION
SM3CP2	\$ 29.00	£ 20.20	€ 25,80	¥ 244.90	SM3 Series Endcap, External Threads

SM3-Compatible Quick-Release Lens Tube Couplers




- Quarter-Turn Connect and Disconnect Adapters
- Lockable Design
- Easily Swap SM3-Threaded Components

NEW
products



Easily Swap Components

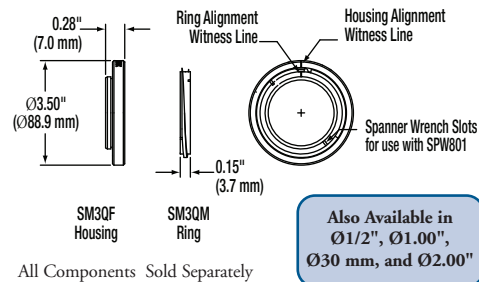

All Components Sold Separately

Also Available in
Ø1/2", Ø1.00",
Ø30 mm, and Ø2.00"

These adapters provide a means to quickly connect and disconnect two similarly SM3-threaded components such as lens tubes, mounted optics, and detectors. By threading a male adapter onto an externally SM3-threaded component and threading the female adapter onto an internally SM3-threaded component, the components can be quickly connected. A witness line on each adapter allows them to be aligned easily. Once they are aligned, a quarter-turn of either component connects the adapters to each other. After the male and female adapters are mated, a nylon-tipped setscrew can be tightened to lock the two adapters together. Similarly, a quarter-turn of either adapter disconnects them. The adapters are available as a set or individually.





Please refer to our website for complete models and drawings.

ITEM#	\$	£	€	RMB	D	DESCRIPTION
SM3Q	\$ 50.40	£ 35.00	€ 44,80	¥ 425.60	Ø3.50" (Ø88.9 mm)	SM3 Quick-Release Adapter Set
SM3QM	\$ 28.00	£ 19.50	€ 24,90	¥ 236.50	—	SM3 Quick-Release Adapter, Male
SM3QF	\$ 28.00	£ 19.50	€ 24,90	¥ 236.50	Ø3.50" (Ø88.9 mm)	SM3 Quick-Release Adapter, Female

SM3 Series Lens Tube Couplers

- Connect to Internally Threaded SM3 Components
- Limit the External Thread Depth of an SM3 Component

Also Available
in Ø1/2",
Ø1.00", and Ø2.00"

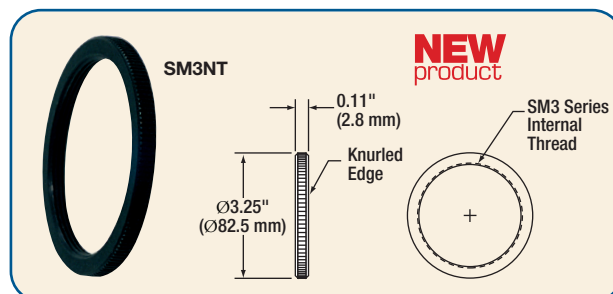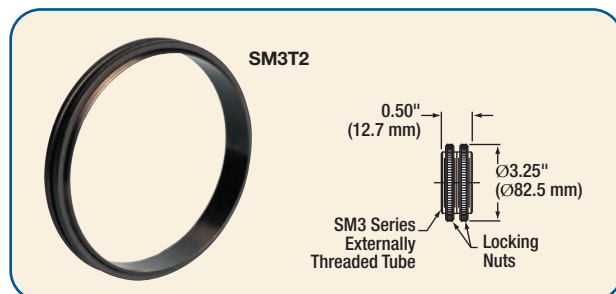

Please refer to our website for complete models and drawings.

ITEM#	\$	£	€	RMB	DESCRIPTION
SM3T2	\$ 37.10	£ 25.80	€ 33,00	¥ 313.30	SM3 Coupler, External Threads
SM3NT	\$ 15.00	£ 10.40	€ 13,40	¥ 126.70	SM3 (Ø3.035"-40) Threaded Locking Nut